

ADVANCED CENTRIFUGES

MAIN ADVANTAGES:

- > Robust Design, QUIET RUNNING
- > Brushless Drive System, MAINTENANCE FREE
- > Automatic Rotor Recognition System, SAFE TO USE
- > Speed Limitation System, SAFE TO USE
- > Microprocessor Controlled Functions of Cover, Time and Speed
- > Adjustable Acceleration and Braking System
- > Special START - STOP Features
- > Fast Spin Features
- > Up to 100 Programs


Centrifuge **CENTRIC 350** is universal laboratory centrifuge, designed for use in medical, industrial and scientific laboratories for separating substances with different relative densities by centrifugal force. The maximum rotational speed of **15.000 RPM** gives the centrifugal force of **21.630 x g**.

A wide usage of programming options and an electronic operation control allow a trouble-free use of the centrifuge. With a special drive, quiet maintenance-free operation without any carbon dust pollution is guaranteed.

Device is equipped with user-friendly options which make the operation and standard settings easier for you. Built-in error detecting functions keep the user from entering incorrect values and check the complete operation.

The centrifuge has the possibility to save programs. In the program store you can save up to 100 different data sets. The centrifuge always keeps the last run program in its memory for an unlimited amount of time allowing the program to be restarted at any time - even if the centrifuge was turned off during operation. All important operation parameters can be seen at a glance.

The settings are executed with knobs and keys on the control panel. The interior of the centrifuge is also easy to clean. We are able to offer you a device that combines functional variety with practical applications.

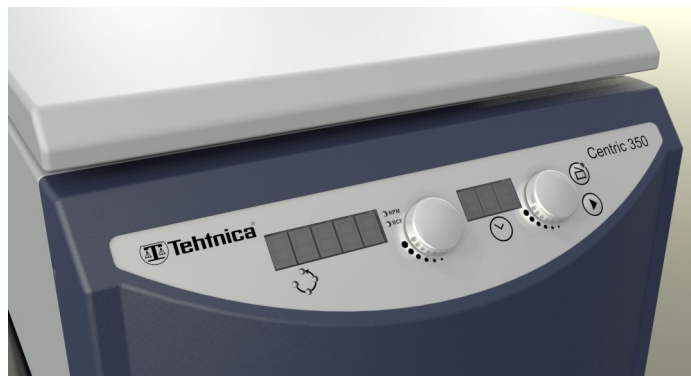

TECHNICAL DATA		CENTRIC 350
Code for	230V 115V	473.000 1473.000
Power supply		230V / 50-60Hz 115V / 50-60Hz
Power requirement		500W
Fuses		2x10A T - 230V 2x16A T - 115V
Isolation category		I
Max. rotation speed		200 to 15.000 rpm, in steps of 10 rpm
Max. centrifugal force		21.630 x g
Max. load		4 x 100ml
Max. kinetic energy		5.020Nm
Max. density of material to be centrifuged		1,2g/ml
Time		0 - 99 min, HOLD function
Program settings		100 programs
Acceleration		0 - 9 levels
Deceleration		0 - 9 levels
Noise level		~ 61dB
Sample temperature at max. speed		+ 15°C
Permitted ambient temperature		+ 2°C to + 35°C
Permitted max. relative air humidity		85%, non condensing
Dimensions (H x W x D)		325 x 405 x 500 mm
Weight		32kg

We reserve the right to alter specification details etc. without prior notice or liability!

Accessories for CENTRIC 350

	FOUR PLACE SWING OUT ROTOR RS 4/100 Code No. 457.200				BUCKET (set of 4) Code No. 462.200			SEALED BUCKET (set of 4) Code No. 464.200			BUCKET (set of 4) Code No. 473.301	
INSERTS												
Insert Code No.	462.203	462.203	462.203	462.203	473.303	462.240	462.208	464.220	462.213	462.218	462.218	462.224
Red. plate Code No. (set of 4)	464.203	464.205	464.207							462.222		
Tube capacity (ml)	0.25	0.5	1.5	5-7	7	5-7	10-12	15	25	50	50	100
Places in insert	12	12	12	12	12	7	6	2	2	1	1	1
Max. No. of tubes	48	48	48	48	48	28	24	8	8	4	4	4
Tube dimension Ø x d (mm)	6 x 30 6 x 45	8 x 32	11 x 40	12 x 75	13 x 105 Vacutainer	12 x 75 Vacutainer	16 x 110 Vacutainer	16 x 120 Falcon	24 x 100	29 x 115 Falcon	36 x 100	44 x 100
Max. tube dimension Ø x d (mm)	6 x 45	8 x 45	11 x 45	12,5x100	13 x 105	12,5x100 cap Ø 17	16,5x105	16,5x121 cap Ø 35	25 x 110	29 x 115 cap Ø 35	36 x 110	44 x 100
Max. radius (mm)	138/151	138	144	151	151	151	148	151	148	148	148	148
Max. speed (RPM)	4400	4400	4400	4400	4400	4400	4400	4400	4400	4400	4400	4400
Max. RCF (g)	2980/3260	2980	3110	3260	3260	3260	3200	3260	3200	3200	3200	3200
Various angle rotors can accept up to 30 tubes arranged in one or two concentric rows. Tubes are placed in angle of 45°. Have in mind that the outer row has greater RCF.	ANGLE ROTOR RA 30/12 Code No. 457.300								ANGLE ROTOR RA 6/50 Code No. 457.320			
												
Shields Code No. (set of 30)	464.280						464.295		464.325 (set of 6)			
Adapter Code No. (1 pcs)	454.141	454.147	454.151			-	-	-				
Adapter capacity (ml)	5	7	10			12	15	50				
No. Of tubes per rotor	30	30	30			30	30	6				
Max. tube dimension Ø x d (mm)	12 x 75 Vacutainer	12 x 100 Vacutainer	16 x 100 Vacutainer			16 x 100 Vacutainer	16 x 120 Falcon	29 x 115 Falcon				
Max. radius (mm)	113 / 69	130 / 115	121 / 106			130 / 115	145 / 130	125				
Max. speed (RPM)	4400	4400	4400			4400	4400	4400				
Max. RCF (g)	2440/1490	2810/2480	2610/2290			2810/2480	3130/2810	2700				
	ANGLE ROTOR RA 24 incl. aluminium lid Code No. 456.600		ANGLE ROTOR RA 30 incl. aluminium lid Code No. 456.620		TWO-PLACE SWING OUT ROTOR Code No. 457.340				MICROTITER PLATE CARRIER (set of 2) Code No. 464.310			
												
Max. capacity (ml)	24 x 2.0		30 x 2.0		No. of plates per rotor				6 plates			
Max. speed (RPM)	15.000		13.000		Size of plate (mm)				86 x 128			
Max. centrifugal force (x g)	21.630		18.510		Max. radius (mm)				132			
Radius (mm)	86		98		Max. speed (RPM)				3000			
Acceleration to max. speed (s)	18		15		Max. RCF (g)				1320			
Braking from max. speed (s)	16		14									
PCR strip ROTOR Code No. 457.380			Vial dimensions Ø x l (mm)		6 x 20							
			Max. no. of vials		48 (6 x 8)							
			Max. speed (RPM)		14.000							
			Max. centrifugal force (x g)		19.500							
			Radius (mm)		86							
			Acceleration to max. speed (s)		15							
		Braking from max. speed (s)		14								

On special request we can provide special custom tube adapter.

* Special dimensioned buckets and inserts

We reserve the right to alter specification details etc. without prior notice or liability!

DOMEL[®]

DOMEL, d.o.o., Otoki 21, 4228 Železniki
BU Laboratory Systems, Na Plavžu 79, 4228 Železniki
 Phone: +386 4 5117 500 Fax: +386 4 5117 501
 WEB: <http://www.domel.com/products/laboratory-equipment>

E-mail: info@tehtnica.si